


Surgiline 550S Surgical Table

Hydraulic

Ref: M2A 000 0 0

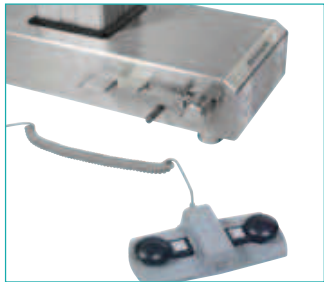

Surgiline 550SE Surgical Table

Electric

Ref: M2B 000 0 0



Bıçakçılar Surgiline 550 series surgical tables are designed to accommodate various type of patient positioning required for demanding surgical procedures with ultimate care to the patient's and user's safety. Surgical Tables are constructed with stainless steel and high quality materials to provide durability and long operating life.

The tabletop consists of four HPL - Compact Laminant sections: head, backrest, seat and leg; each with removable, easy to clean, anti-bacterial, and waterproof polyurethane pads for patient comfort. The radiographic tabletop* made of HPL - Compact Laminant can be easily installed and removed for X-ray procedures.

The durable cast iron table base is covered with a stainless steel. The "T" shape design of the Surgiline 550 series allows easy access to the tabletop in any position. Special brakes and anti-static casters allow ease of mobility. Stainless steel side rails accept a variety of available accessories.

Surgiline 550 series surgical tables are reliable products having both quality and competitive pricing it is an excellent choice.

Surgiline 550S Surgical Table, Hydraulic

The tabletop can be raised and lowered to the desired height with a hydraulic lift in a quiet and non-vibrating manner by using a foot pedal. Downward speed is adjustable. Other specific table settings are performed manually. The tabletop can be rotated for 360°.

Surgiline 550SE Surgical Table, Electric

The tabletop can be raised and lowered to the desired height with an electrically operated piston in a quiet and non-vibrating manner. This is achieved with a foot pedal control. Other specific table settings are performed manually.

(*) *Optional*

Technical Details

Height adjustment:

770 - 1000 mm (550S)
815 - 1115 mm (550SE)

Height adjustment (w/mattress):

800 - 1030 mm (550S)
845 - 1145 mm (550SE)

Tabletop (w/mattress):

465 x 1950 mm

Tabletop (w/o mattress):

755 x 1950 mm

Table base: 550 x 800 mm

Head section: 465 x 255 mm

Backrest section: 465 x 550 mm

Seat section: 465 x 445 mm

Right leg section: 230 x 645 mm

Left leg section: 230 x 645 mm

Electro-mechanical piston

(550SE): DC 24V

Packing Details

Net weight:

220 kg (550S)

180 kg (550SE)

Gross weight:

302 kg (550S)

262 kg (550SE)

Packing dimensions:

1000 x 1680 x 1100 mm

Positions

Head section: Upward 20°, downward 90° Extendable 5 cm, Removable

Lateral: Right 20°, left 20°

Trendelenburg: 25°

Reverse Trendelenburg: 25°

Backrest section: Downward 15°, Upward 75° (w/mattress), Upward 90° (w/o mattress), Upward 20° (w/radiographic tabletop and mattress)

Reflex: 105°

Flex: 195°

Leg section: Downward 90°, Upward 15° Spreading 90°, Removable, Independent movement

Renal surgery: 165°

Tabletop: Rotates 360° (optional in 550SE)

