

BA8800 PamaBed



The
Excellent Guide

Towards Patient Care



Hospital ICU/CCU Electrical Bed

www.lklbeds.com



TUV Product Service

The model BA8800 series was introduced as an excellent guide towards patient care. Fully succeeding the quality standards of the ISO 13485:2003, it is fully electric-motorized thus making it user-friendly. A combination of innovation and state-of-art technology the BA8800 is all about patient comfort redefined and it is a complete solution for caregivers.

> X-Ray Translucent Base

X-Ray Translucent base from head to umbilicus completed with x-ray cassette holder

> Patient Control embedded Side Rail

To allow patient to have desired & comfortable position

BA8

> Side Rail

Split version PU side rails adjustment with gas spring

> Tente Twinwheel


High quality 150mm dia. Tente Twinwheel with central locking system

> Drainage Bag Holder

> Manual override CPR Release of both side for backrest

> ABS Mattress Base

Mattress base tops made of ABS (high impact plastic) removable, detachable and easy to clean

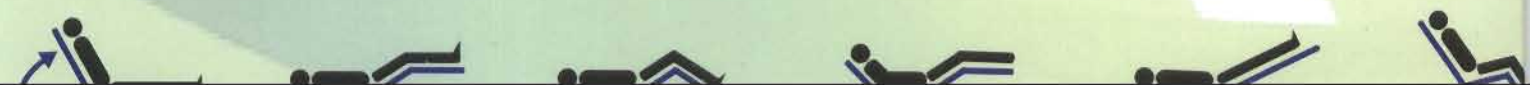

800 PamaBed



> Nurse Lock

Caregiver control

- Nurse lock control located at foot end of the bed. Featuring with CPR one press button for neutral position

:: Features

- + All parts made by tooling to ensure uniform quality
- + Steel bed frame with scratch and chip resistant epoxy coating
- + PE Head & Foot boards with locking mechanism for easy attachment, detachment and can be cleaned easily
- + Split fold down PE Integral safety side rails with convenient one-hand operation for engaging and disengaging
- + X-Ray radio Translucent base from head to umbilicus section complete with X-Ray cassette holder tray with adjustable film positions
- + Mattress base tops made of ABS (high impact plastic) easy removable and detachable
- + Four high-quality 150mm diameter Tente Twinwheel with central locking, free and steer system with pedal located at foot end of the bed for convenient operation.
- + Hi-Lo, Backrest, Legrest, Trendelenburg and Reverse Trendelenburg by electrical motorized
- + Embedded Patient Controls button on inner side rails of the head section to allow patient to have own desired and comfortable position
- + Caregiver control (nurse lock and nurse control) located at foot end of the bed
- + Nurse lock to lock out the individual function of the side rails
- + CPR one press button (on handset controls located at foot end of the bed) for CPR control to neutral the bed position (lowest and horizontal position)
- + Manual Override Quick CPR Release Lever for lowering the backrest located at both sides of the bed
- + Auto contour function allows patient to achieve the most ergonomic and comfortable 'cardiac chair position', while reducing the need to move patient to alternative seating
- + Provisions of 6 IV Pole attachment and 2 monkey pole attachment
- + Drainage bag holder located at both sides of the bed
- + PU bumper at 4 corners
- + IV Drip Pole storage

:: Electrical System

- + Linak DC-actuator system (TUV & UL approval & EC60601-1)
- + 24 VDC with built in back-up battery
- + Input Voltage: 220 - 240V; 50Hz - 60Hz
- + Safety Rating
Controllers: IP54, Actuators: IP54, Nurse Control: IP66

:: Technical Specifications

- + Backrest: 0° - 80°
- + Kneerest: 0° - 35°
- + Height Adjustment(mm): 480-780 (floor to mattress base)
330 (From mattress base to PE Side Rail)
- + Trendelenburg & Reverse Trendelenburg: -15° to +15°
- + Mattress Dim(mm): 1980(L) x 915(W) x 100(H)
- + Overall Dim(mm): 2200(L) x 1050(W) x 480 - 780(H)
- + Safety load capacity: 230kg UDL

:: Accessories Included

- + Hospital Medi-Flex Mattress with memory foam to allow patients' body contouring and reduce pressure, mattress cover of fire retardant, anti-bacterial and fluid resistant (x 1 unit)
- + Stainless steel Telescopic IV Drip Pole with 2 hooks (x 1 unit)

:: Optional Accessories

- + Monkey pole
- + Oxygen cylinder holder





278577-X

LKL Advance Metaltech Sdn. Bhd.
 No.3, Jalan BS 7/18,
 Kawasan Perindustrian Bukit Serdang,
 Seksyen 7, 43300 Seri Kembangan,
 Selangor Darul Ehsan,
 Malaysia.
 Tel: +603 8948 2990 (Hunting Line)
 Fax: +603 8948 7904, 8943 8487
 Email: info@lklbeds.com
 URL: www.lklbeds.com

Achievements & Awards :

