


BA9000



The
Excellent Guide

Towards Patient Care

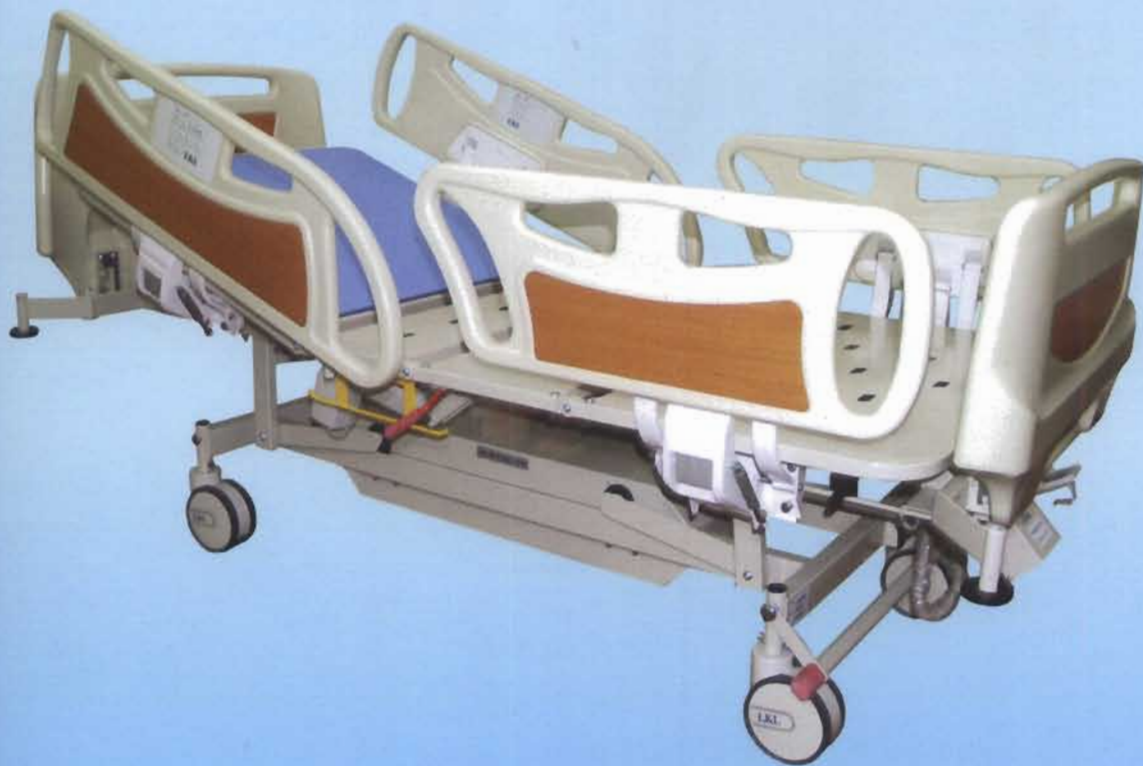

www.lklbeds.com

Hospital ICU / CCU Electrical Bed



TUV Product Group

BA9

> Side Rail

Integral safety side rail of split version
- 4 individual PU side rails with gas spring assisted adjustment

> X-Ray Translucent Base

X-Ray Translucent base from head to umbilicus complete with x-ray cassette holder

> Patient Control Embedded Side Rail

To allow patient to have desired & comfortable position

> Manual Override CPR release lever of both side for backrest

> Drainage Bag Holder

> ABS Mattress Base


Mattress base tops made of ABS (high impact plastic) removable, detachable and easy to clean

> Tente Twinwheel

High quality 150mm dia. Tente Twin-wheel with central locking system

> Nurse Lock

Caregiver control
- Nurse lock control located at the end of the bed. Featuring one press button to the desired position or back to neutral

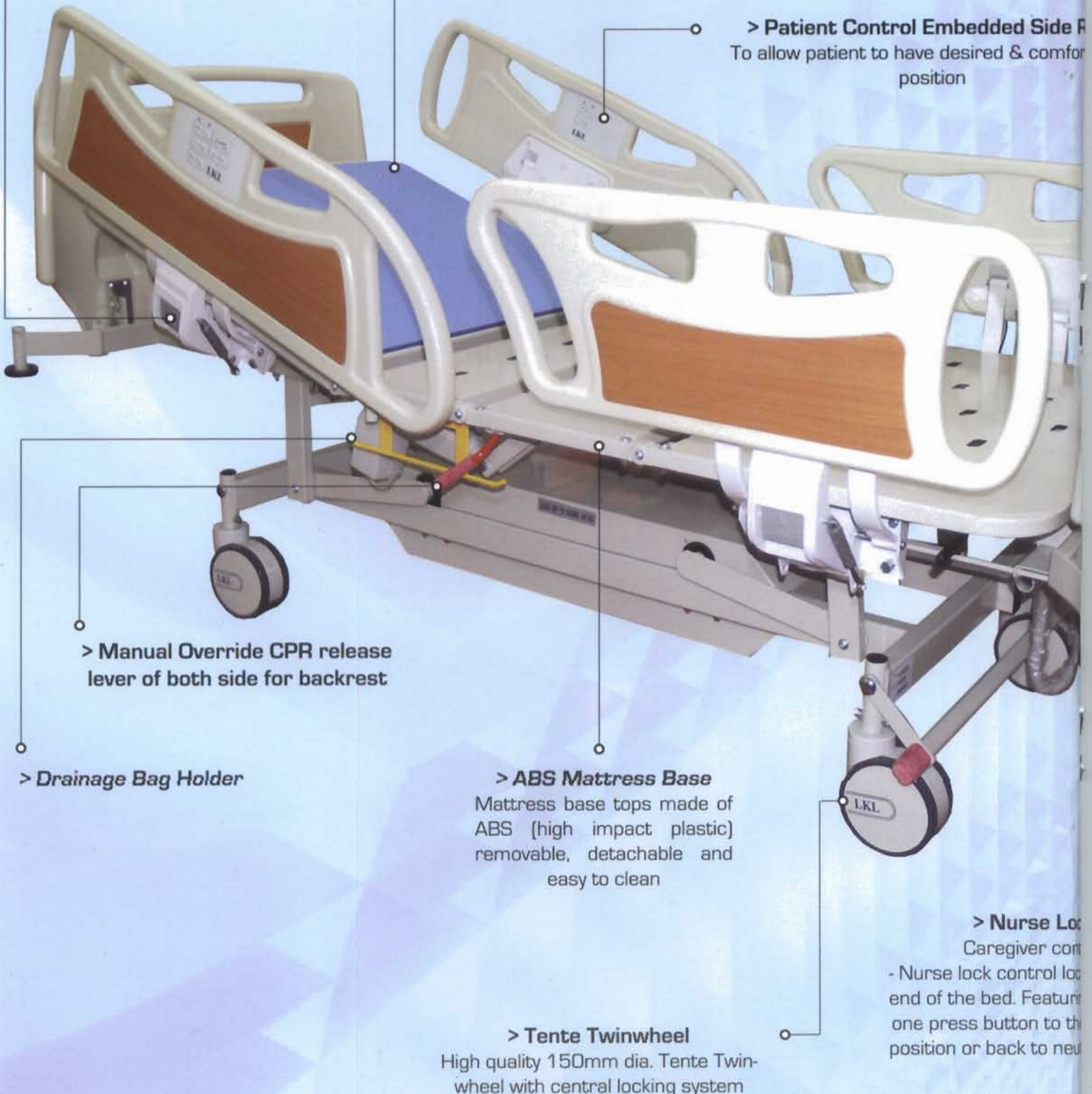

000

il
ble



ed at foot
with CPR
horizontal
al position

:: Features

- + All parts made by tooling to ensure uniform quality
- + Steel bed frame with scratch and chip resistant epoxy coating
- + PE Head & Foot boards with locking mechanism for easy attachment, detachment and can be cleaned easily
- + Split fold down PE Integral safety side rails with convenient one-hand operation for engaging and disengaging
- + X-Ray radio Translucent base from head to umbilicus section complete with X-Ray cassette holder tray with adjustable film positions
- + Mattress base tops made of ABS (high impact plastic) easy removable and detachable
- + Four high-quality 150mm diameter Tente Twinwheel with central locking, free and steer system with pedal located at foot end of the bed for convenient operation.
- + Hi-Lo, Backrest, Legrest, Trendelenburg and Reverse Trendelenburg by electrical motorized
- + Embedded dual sided Patient Controls button on both outer and inner side rails of the head section
- + Full function on the outer side (Nurse Controls button) for caregiver and limited function on the inner side (Patient Controls button) to allow patient to have own desired and comfortable positions
- + Caregiver control (nurse lock and nurse control) located at foot end of the bed
- + Nurse lock to lock out the individual function of the side rails
- + CPR one press button (on handset controls located at foot end of the bed) for CPR control to neutral the bed position (lowest and horizontal position)
- + Manual Override Quick CPR Release Lever for lowering the backrest located at both sides of the bed
- + Auto contour function allows patient to achieve the most ergonomic and comfortable 'cardiac chair position', while reducing the need to move patient to alternative seating
- + Provisions of 6 IV Pole attachment and 2 monkey pole attachment
- + Drainage bag holder located at both sides of the bed
- + PU bumper at 4 corners
- + IV Drip Pole storage

:: Electrical System

- + Linak DC-actuator system (TUV & UL approval & EC60601-1)
- + 24 VDC with built in back-up battery
- + Input Voltage: 220 - 240V; 50Hz - 60Hz
- + Safety Rating
Controllers: IP54, Actuators: IP54, Nurse Control: IP66

:: Technical Specifications

- + Backrest: 0° - 80°
- + Kneerest: 0° - 35°
- + Height Adjustment(mm): 480-780 (floor to mattress base)
360 (From mattress base to PE Side Rail)
- + Trendelenburg & Reverse Trendelenburg: -15° to +15°
- + Mattress Dim(mm): 1980(L) x 915(W) x 100(H)
- + Overall Dim(mm): 2200(L) x 1050(W) x 480 - 780(H)
- + Safety load capacity: 230kg UDL

:: Accessories Included

- + PatientMate Technology Mattress with memory foam to allow patients' body contouring and reduce pressure, mattress cover of fire retardant, anti-bacterial and fluid resistant [x 1 unit]
- + Stainless steel Telescopic IV Drip Pole with 2 hooks [x 1 unit]

:: Optional Accessories

- + Monkey pole
- + Oxygen cylinder holder





278577-X

LKL Advance Metaltech Sdn. Bhd.
No.3, Jalan BS 7/18,
Kawasan Perindustrian Bukit Serdang,
Seksyen 7, 43300 Seri Kembangan,
Selangor Darul Ehsan.
Malaysia.

Tel: +603 8948 2990 (Hunting Line)

Fax: +603 8948 7904, 8943 8487

Email: info@klbeds.com

URL: www.klbeds.com

Achievements & Awards :



MPC Productivity Award
2007



Woman Entrepreneur
2007



MMEBE Platinum
Brand Award



MMEBE Corporate Branding
Export Market (Gold)



5th Malaysia
100 Outstanding SMEs



50th Enterprise Award
Ranked No.6



SMB Emerging Business



SMI - Maybank Rising Star