

Hospital Hydraulic Bed

Hospital Hydraulic Bed

The *Excellent Guide*

Towards Patient Care



www.lklbeds.com



BA 4120-T



Foot Pedal
Controls HiLo adjustable hydraulic jack

BA4120-T Double Fowler

Product Features

- + Steel construction
- + HiLo adjustable by hydraulic jack controlled by foot pedal
- + Backrest & legrest adjustable by gas spring mechanism
- + Trendelenburg/reverse trendelenburg assisted by gas spring mechanism
- + Detachable ABS head and foot boards
- + 150mm central locking wheels with steer, free & total brakes function
- + Inclusive of collapsible aluminium siderail and medi-flex mattress
- + 2 provision for IV Drip Pole attachment
- + Safety load capacity: 230kg



Central Locking System
150mm castors with Steer, Free & Total Brakes

Bed Positions

+ Backrest (0° - 80°) and
Legrest (0° - 40°)

+ Hi-Lo Position from Floor
to Mattress Support Base
(460 mm - 820 mm)

+ Trendelenburg (0° - 12°) and
Reverse Trendelenburg (0° - 6°)




Hospital Hydraulic Bed

Ηοςπιτάι Ηλδραυλικ Βεδ


The BA 4000 series was developed with the healthcare personnel and patients in mind. It provides a safer and smarter alternative to regular hospital beds. Its advanced hydraulic and gas spring system allows care givers to operate the bed with ease, hence eliminating occupational hazards and enhancing the overall quality of patient care.

The BA 4000 series... a smarter choice...

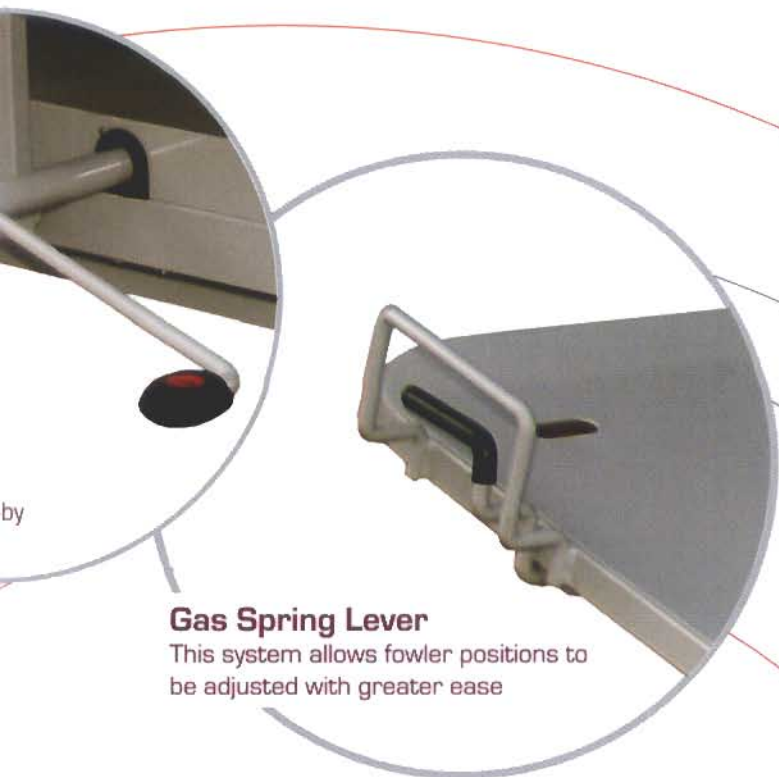

Gas Spring Lever

This system allows fowler positions to be adjusted with greater ease

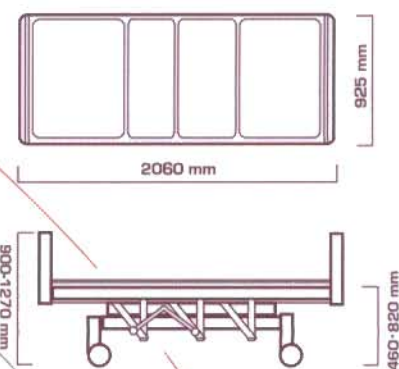

Features

tem
r, Free & Brake system



LKL Beds utilize only the highest quality "enabling technology" components. Among the relevant standards and



VDA






BA 4121-T Double Fowler

Features

- + Steel construction
- + Hilo adjustable by hydraulic jack controlled by foot pedal
- + Backrest adjustable by gas spring mechanism & retractable crank legrest
- + Trendelenburg/reverse trendelenburg assisted by gas spring
- + Detachable ABS head and foot boards
- + 150mm central locking wheels with steer, free & total brakes function
- + Inclusive of collapsible aluminium siderail and medi-flex mattress
- + 2 provision for IV Drip Pole attachment
- + Safety load capacity: 230kg

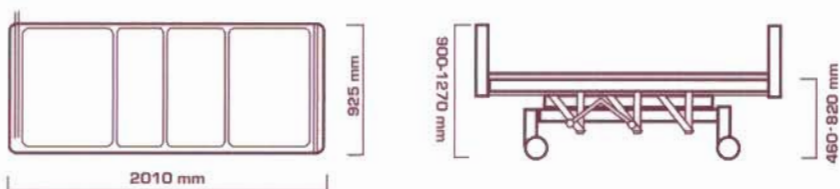

BA4000 Series Beds

	Head/Foot Board	Fowler Position		Fowler Mech nism				Special Features	
	ABS	Single	Double	Head		Leg		Trendelenburg	X-Ray Translucent Base
				Crank	Gas Spring	Crank	Gas Spring		
BA 4110	●	●			●				
BA 4120	●		●		●		●		
BA 4121	●		●		●	●			
BA 4120-T	●		●		●		●	●	
BA 4121-T	●		●		●	●		●	
BA 4120-TX	●		●		●		●	●	●
BA 4122-T	●		●		●	●		●	

Product Features




Bed Stripper

For Linen or personal items placement



Bed Tilting Levers

For Linen or personal items placement

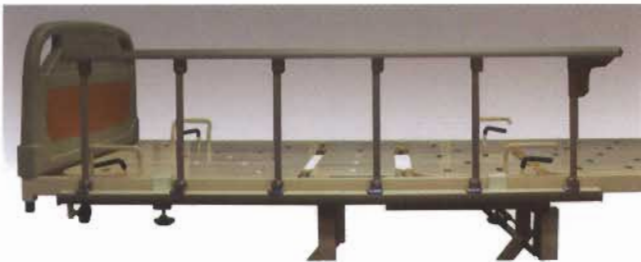

BX200

Collapsible Aluminium Siderail

- + ABS top
- + Dim(mm): 410(H) to 1440(L)



BX411

Medi-Flex Mattress

Mattress foam:

- + PU foam with 32 density

Mattress cover:

- + Water-Proof
- + Fire Retardant
- + Anti-Bacterial Treated
- + Anti-Static
- + Stain Retardant
- + Urine Resistant

- + Dim(mm): 1905(L) x 889(W) x 102(H)



Central Locking Castors

150mm diameter Tente (Germany) castors with unrivalled durability and mobility.



278577-X

LKL Advance Metaltech Sdn. Bhd.
No.3, Jalan BS 7/18,
Kawasan Perindustrian Bukit Serdang,
Seksyen 7, 43300 Seri Kembangan,
Selangor Darul Ehsan,
Malaysia.

Tel: +603 8948 2990 (Hunting Line)

Fax: +603 8948 7904, 8943 8487

Email: info@klbeds.com

URL: www.klbeds.com

Achievements & Awards :



MPC Productivity Award
2007



Woman Entrepreneur
2007



MMEBE Platinum
Brand Award



MMEBE Corporate Branding
Export Market (Gold)



5th Malaysia
100 Outstanding SMEs



50th Enterprise Award
Ranked No.6



SMB Emerging Business



SMI - Maybank Rising Star