


Surgical Blades and Handles Range



Swann-Morton®

for both up and down cutting. Individually packed, sterilised and boxed in 10's. This blade fits the 'FINE' range of handles.

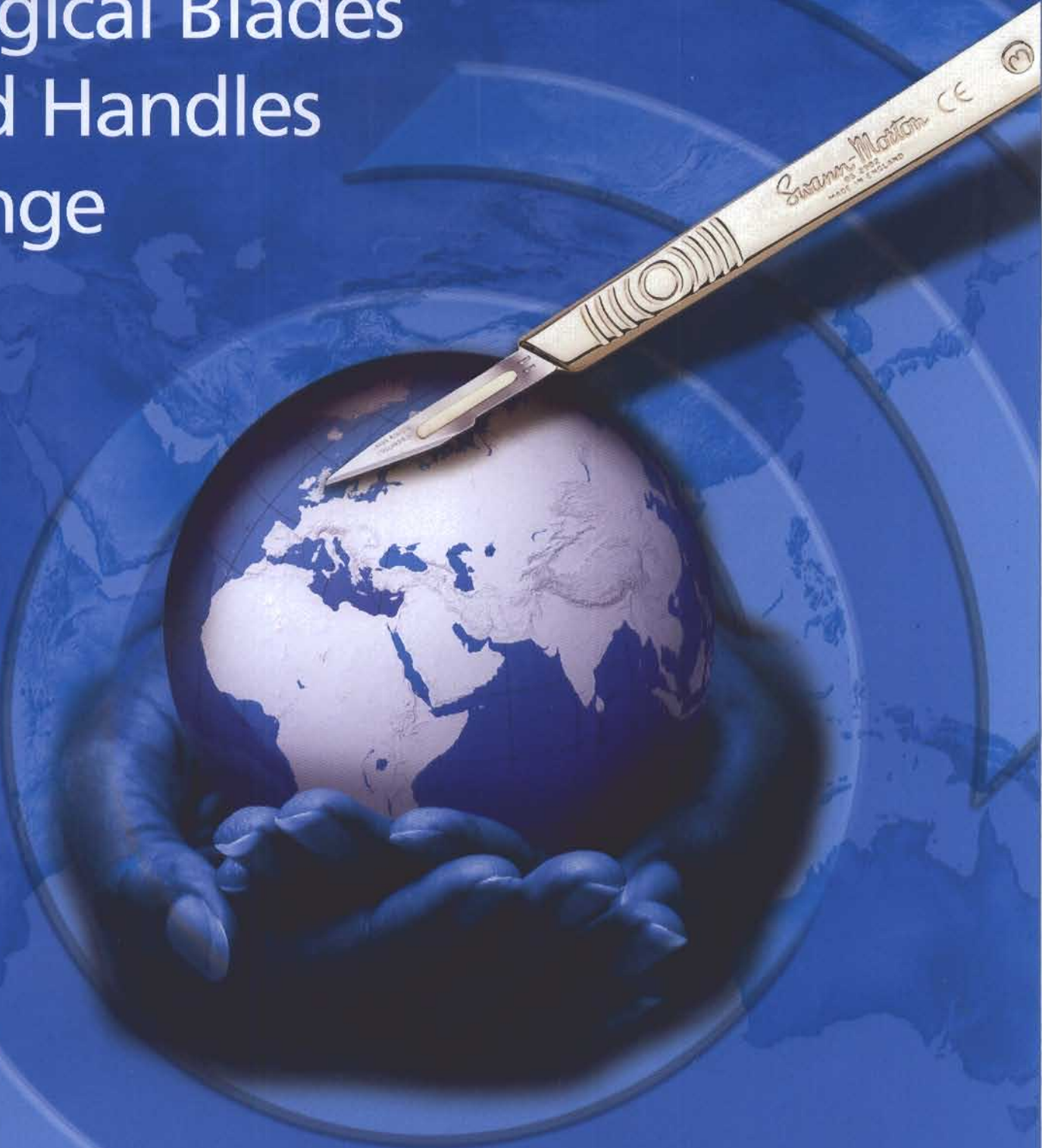

Surgical Blades and Handles by

Swann-Morton®



PM40 (Stainless Steel)

PM8



PM40B

PM60





The Quality Edge

Swann-Morton scalpels and blades have a reputation for quality worldwide that is second-to-none.

That reputation has been built on exceptionally – and consistently – high standards of durability and sharpness. Swann-Morton is accredited to BS EN ISO 13485, enabling the company to CE mark its products under the Medical Device Directive (93/42/EEC). All standard fitting surgical blades and handles are produced in accordance with BS 2982, ISO 7740 / BS EN 27740. Every single batch of blades is subject to the strictest possible quality control procedures at every single stage of its manufacture, and the company also has a fully equipped microbiological laboratory.

Such attention to quality is, of course, vital for Swann-Morton. The company's products play an important role in one of the most demanding environments in the world – the operating theatre.

For the same reason, Swann-Morton blades are subject to constant development and refinement. As surgical techniques become ever more highly refined, so must Swann-Morton develop blades that are equal to the challenge.

Surgeons in just about every country in the world trust other people's lives to Swann-Morton blades. There can surely be no better evidence of the importance the company places on quality.

for both up and down cutting. Individually packed, sterilised and boxed in 10's. This blade fits the 'FINE' range of handles.

Surgical Blades

Swann



6



4L



Graduated
4



4



L (Long
Handle)



S (Standard
Handle)



Cervical Biopsy Blade

In stainless steel, sterilised, supplied in single plastic trays in boxes of 10. Fits both Swann-Major handles.



PM40

PM40B



PM60

and Handles by

Morton



PM40
(Stainless Steel)

PM8

Standard surgical blades in carbon steel (sterile or non-sterile) or in stainless steel (sterile). Boxed in 100's (sterile blades individually foil-wrapped; non-sterile sealed in foil packs of 5). Non-sterile blades also available singly packed in peel-open packets, boxed in 100's. Sabre E/11 and D/15 blades are produced in carbon steel only, sterile, boxed in 100's.



'Fine' Blades and Handles

The 'Fine' range has small, slim, precise blades in stainless steel for delicate operations. Six special handles, including the SF13 and SF23 with independently screwed top sections to facilitate blade-changing, take all blades, supplied in individual metal foil packets, available in dispenser units of 25. Handles are supplied individually in plastic wallets.



Swann-Major Surgical Blades and Handles

All these extra-strong blades fit either handle. Blades are carbon steel, sterile, individually foil-wrapped, boxed in 50's.



Post-Mortem Knife - PM40 and PM40B

Post-mortem knife with a choice of heavy-duty carbon steel blades. Screw-type fitting to stainless steel handle. Spare blades, individually wrapped, boxed in 10's.



Post-Mortem PM8 Handle & PM60 and PM60B Blades

Complementing existing Post-Mortem range, offering easy blade change, in non-sterile carbon steel. Both blades VPI protected in single-peel packets, boxed in 50's. Handles supplied singly in plastic wallets.



Stitch Cutters

Requiring no handles, the standard stitch cutter is in carbon steel, the midi and long versions in stainless steel.

Available individually foil-wrapped in sterile packs of 100.

Sterilisation is by Gamma Radiation (minimum dosage 25kGy, 2.5 Mrad) carried out under strict laboratory control in Swann-Morton (Services) Ltd's Cobalt 60 continuous irradiation plant.

All items are shown actual size



All medical devices are CE marked in accordance with the Medical Device Directive (93/42/EEC).



FM 73368

"Swann-Morton" and the 'ring pattern' logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustrations used in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in



The Swann-Morton Surgical Blade Remover

Swann-Morton's Surgical Blade Remover is for single use, and is available sterile in boxes of 50, and non-sterile in boxes of 100. It is designed to remove surgical blades from standard No. 3 and 4 fitment handles, and protects hospital staff removing scalpel blades after operations, guarding against cuts and abrasions.



Skin Graft Blade (SG3)

In stainless steel, with 5cm long cutting edge. Fits a No.3 fitment handle. Sterilised, packed individually, boxed in 20's.



Stitch Cutter

Stainless steel with No.3 handle fitment. Individually foil-wrapped in boxes of 100.





16
40

B3L

B3

9

7

5

3L

Morton CE

Swann-Morton
CE 5

Swann-Morton

Myringotomy

'Minor' Disposables

Can be used as part of a "Ward Kit". Small plastic handle can be supplied fitted with blade No's 10, 11, 15 or 15A. Polythene blade guard. Non-sterile and bulk-packed.



Sabre E/11



Disposable Scalpels

Available with all standard stainless steel blade shapes, complete with plastic handles, individually wrapped in sterile packs of 10.



Guarded Disposable Scalpels

With an ultrasonically-welded blade fitment, these scalpels are available with all standard stainless steel blade shapes, and are fitted with a polythene protective guard. The sterile scalpels are individually wrapped in a view pack, and boxed in units of 10.

Also available in a bulk, non-sterile format.



Sabre Disposable Scalpels

A choice of 3 'Sabre' sterile disposable scalpels – E/11, D/15 and B/23 – each supplied individually packed, boxed in 10's.



Disposable Scalpels with Retractable Blade

These disposables have been specially designed to reduce accident risk potential by having a high-visibility red handle, retractable blade operated by a colour-coded button which audibly clicks to indicate open and locked position, and ribbed, slip-resistant finger grips for secure hold in gloved hands.

They are available with stainless steel blade No's 10 (black button), 11P (blue button) and 15 (green button), in sterile packs of 25.


Also available in non-sterile packs of 500.



Standard Surgical Blades and Handles – No.3 Fitment

Surgical Blade No's 6, 9, 10, 10A, 11, 11P, E/11, Sabre E/11, 12, 12D, 13, 14, 15, 15A, 15C, 15T, Sabre D/15, 16, 40 fit No's 3, 3L, 5, 7, 9, B3 & B3L handles.

Handles are available in nickel, with selected sizes available in stainless steel. All handles are individually wrapped. The easily gripped B3 & B3L are light, adaptable stainless handles, supplied individually in plastic wallets.

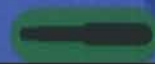

Standard Surgical Blades and Handles – No.4 Fitment



For further information contact:

Swann-Morton

Swann-Morton Limited,
Penn Works, Owlerton Green, Sheffield S6 2BJ.
Tel: 0114 234 4231/234 4223 (Sales)
Fax: 0114 231 4966

e-mail:
uksales@swann-morton.com
exportsales@swann-morton.com

